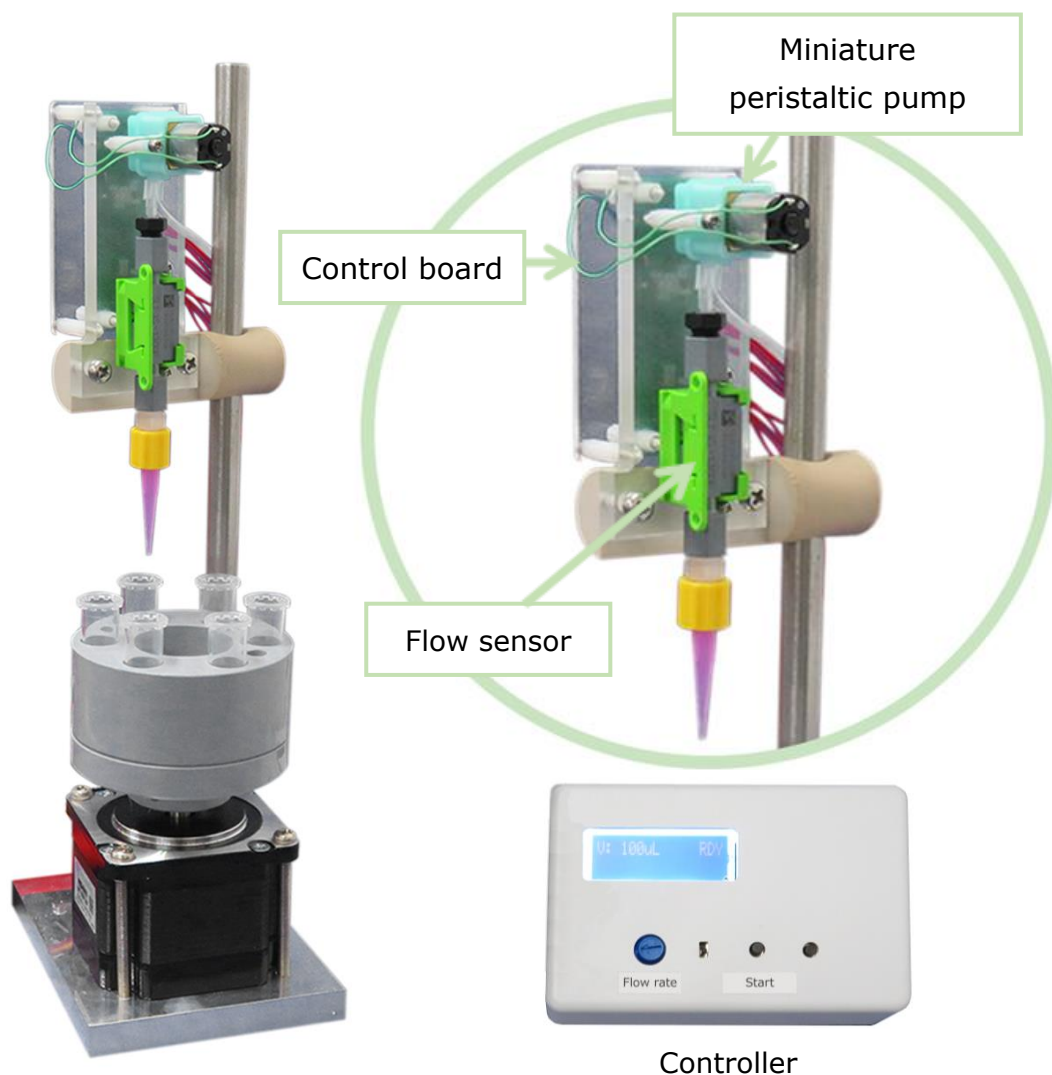


INTEGRATED PRODUCT

Smart & Easy Dispensing Module**High-precision dispensing at low cost****Features**

- High-precision dispensing can be carried out easily even with a peristaltic pump, which has a pulsating flow. Simply enter the target dispensing volume into the controller, and the high-precision flow sensor monitors the total flow to give feedback-control to the pump. No need to control an expensive and complicated syringe pump anymore.
- The control unit can be used in any combination with the various types of pumps that we offer. As an example, the ultra-small model shown in the picture above can achieve high precision of $\pm 1\%$ with a dispensing volume of 300 μl or more.

TFS (Takasago Fluidic Systems)

1900 West Park Drive, Suite 280, Westborough, MA 01581 U.S.A.

Tel: 508-983-1434 FAX: 508-983-1401

E-mail: info@takasago-fluidics.com URL: <http://www.takasago-fluidics.com/>